



DEPPELER

Signature of excellence

SHARPENING
SOLUTIONS




WHY CHOOSE INSTRUMENTS THAT CAN BE RESHARPENED?

A study tested and compared Deppeler instruments and instruments that cannot be resharpened. The results show that, after 1000 strokes on dentin and five reprocessing cycles, all instruments have lost about 75% of their effectiveness.

At this stage, if the practitioner uses an instrument that cannot be resharpened, the blunt instrument will continue losing its efficiency.

If the practitioner uses a Deppeler instrument, the cutting edge can be **resharpened to regain optimal sharpness and comfort.**

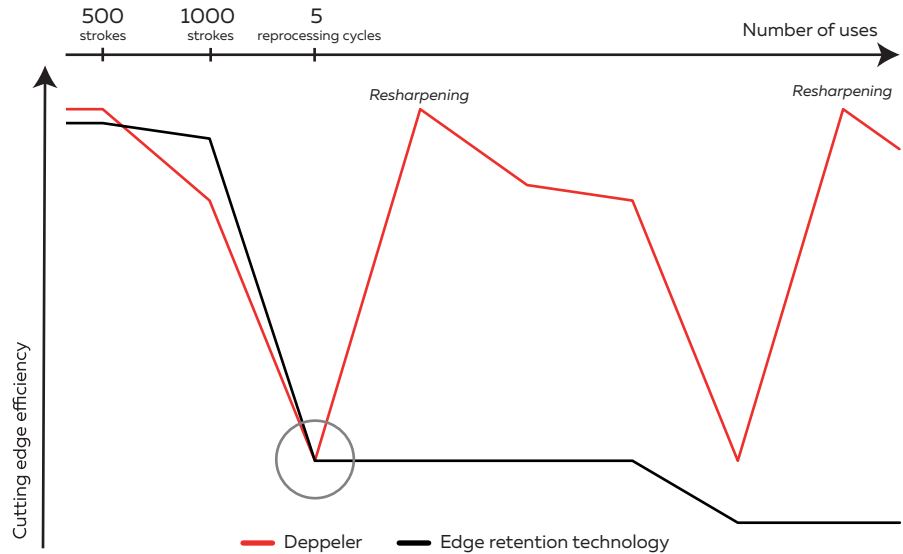

Schmidlin et. al.: In vitro evaluation of 3 curettes with edge retention technology after extended use. SSO Journal 119: 1200-1207. (12/2009)

Flexibility

The working parts of resharpenable instruments are more flexible and less likely to break in the patient mouth.

Efficiency

A bare steel edge is more efficient than a coated edge.

Savings

Resharpenable instruments need to be replaced less frequently. Their cost is lower in the long run.

Sustainability

In addition to having to purchase instruments less frequently, all-steel instruments are also easier to recycle.

SHARPENING METHODS

Maintain the cutting edges of your scalers and curettes until the end of their life cycle. With Deppeler's foolproof Easy Sharp solutions, sharpen your instruments easily and effectively, without any facets, serrated edges, or burrs.




MANUALLY

Natural stones



Arkansas Hard Black
ESS1
For regular sharpening

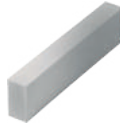

Arkansas Translucent
ESS3A
Finest quality
For faster sharpening
and better results



Silicone Carbide
ESS4
For titanium instruments
and very blunt steel
instruments



WITH A COACH

Easy Sharp Coach
ESC

Designed for visual guidance and easy control of the sharpening angle.
1 coach for universal curettes and scalers, 1 coach specifically for Gracey curettes.
Suitable for all brands.

How to use an
Easy Sharp Coach?



WITH A MACHINE

Easy Sharp Device
ESD

Designed for maximum precision of the sharpening angle and optimal control of the movement.
The cutting edge feels brand new, back to its original state.
Suitable for Deppeler instruments.

How to use an
Easy Sharp Device?



Order information

EASY SHARP COACH

ref. ESC

includes:

- 2 coaches (1 universal and 1 Gracey)
- 1 Arkansas Hard Black stone ESS1

ref. ESCPLUS

includes:

- 2 coaches (1 universal and 1 Gracey)
- 1 Arkansas translucent stone ESS3A

EASY SHARP DEVICE

ref. ESD

includes:

- Device
- 1 Arkansas Hard Black stone ESS1

ref. ESDPLUS

includes:

- Device
- 1 Arkansas translucent stone ESS3A

We recommend resharpener your instruments
once every 5 to 8 treatments.



WWW.DEPELER.CH

f DEPPELERSA | @ DEPPELER_INSTRUMENTS



DEPELER

Signature of excellence